## Personenwaage Victoria

Victoria Bathroom Scales | Báscula personal Victoria

Pèse-personne Victoria | Bilancia pesapersona Victoria

Personenweegschaal Victoria | Waga osobowa Victoria

Deutsch

English

Español

Français

Italiano

Nederlands

Polski

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BM 701 + 708

Operating instructions | Instrucciones de operación | Mode d'emploi

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# Dear Customer!

You have chosen to purchase a high-quality product from the brand ADE, which combines intelligent functions with an exceptional design. These bathroom scales make it possible for you to determine and monitor your body weight. The many years of experience of the brand ADE ensures a high technical standard and proven quality.

We hope you enjoy weighing yourself and wish you success in reaching your goal weight!

Your ADE Team

## General Information

---

### About this manual

This operating manual describes how to safely use and care for your products. Keep this operating manual in a safe place in case you would like to refer to it later. If you pass this item on to someone else, pass on this manual as well.

Failure to observe the recommendations in this operating manual can lead to injuries or damage to your product.

### Explanation of symbols

This symbol, combined with the word DANGER, warns against serious injury.

This symbol, combined with the word WARNING, warns against moderate and mild injuries.

This symbol, combined with the word NOTE, warns against material damage.

This symbol signifies additional information and general advice.

# Table of Contents

---

General Information .....	8
Intended Use.....	10
Safety .....	10
At a glance (scope of delivery) .....	11
Calibrating .....	12
Weighing .....	12
Cleaning.....	12
Fault/Remedy.....	12
Technical Data .....	13
Guarantee .....	13
Disposal .....	13

## Intended Use

---

The bathroom scales are designed for the weighing of people in the home. It is not suitable for use in clinics, medical practices, etc.

The maximum carrying capacity of 136 kg must not be exceeded.

## Safety

---

- These scales can be used by children 8 years of age and above, as well as by persons with reduced physical, sensory, or mental capacities or lacking experience and knowledge, if they are supervised or have been instructed in the safe use of the scales and understand the dangers that may arise.
- Children must not play with the scales.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.

### **DANGER for children**

- Keep children away from packaging material. Risk of choking if swallowed. Also beware that children do not pull the packaging material over their heads.

### **WARNING against risk of injury**

- Do not step onto the outer edge of any one side of the scales. Risk of tipping!
- Do not step onto the scales with wet feet or when its surface is damp. Risk of slipping!


### **BEWARE of material damage**

- Place the scales on a stable, even surface.
- Do not overload the scales. Pay attention to the "Technical Data". Place yourself on the scales exclusively for weighing purposes.
- Protect the scales from sharp jolts and vibrations.
- Do not make any modifications to the scales. Only allow repairs to be carried out at a specialist workshop. Repairs carried out incorrectly can considerably endanger the user.

- There are feet on the underside of the scales. Since floor surfaces are made from various materials and are treated with diverse care products, the possibility cannot be completely excluded that some of these substances contain constituents that attack and soften the feet of the device. If necessary, place an anti-slip underlay beneath the scales.

## At a glance (scope of delivery)

---

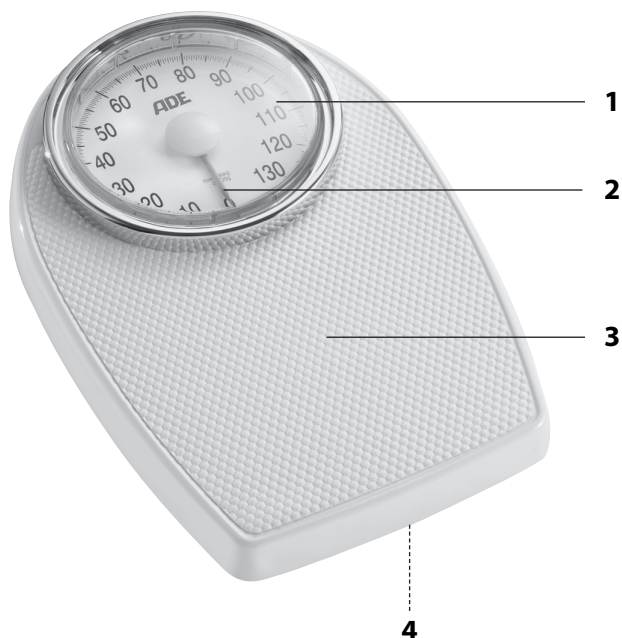- 1 Dial
- 2 Weight pointer
- 3 Weighing surface
- 4 Calibration wheel

Not depicted: Operating instructions

# Calibrating

---

Before weighing yourself for the first time, you have to calibrate the scales in order to obtain a weighing result that is as exact as possible.

- Using the calibration wheel, centre the weight pointer on the “0” of the dial.

Check the zero-setting from time to time, e.g. when you change the location of the scales. When no load is on the scales, the weight pointer must be centred on “0”.

# Weighing

---

1. Position yourself on the scales. Remain standing still and do not hold on to anything.
2. Read off your measured body weight from the weight pointer.

- Weigh yourself with as little clothing on as possible.
- Try to weigh yourself always at the same time and under consistent conditions in order to obtain comparable values.

# Cleaning

---

## NOTE on material damage

- Do not immerse the scales in water or other liquids.
  - Do not under any circumstances use harsh, scraping, or abrasive cleaning agents or the like to clean the scales. If you do, the surface could be scratched.
- Wipe down the scales as needed using a lightly dampened, soft cloth.
  - Disinfect the weighing surface from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

# Fault/Remedy

---

Fault	Possible Cause
The scales are showing obviously false weighing results.	<ul style="list-style-type: none"><li>▪ Are the scales standing on a firm, level surface?</li><li>▪ Are the scales correctly calibrated? (see chapter “Calibrating”).</li></ul>

## Technical Data

---

Model:	Victoria Bathroom Scales BM 701 (white) BM 708 (black)
Measuring range:	max. 136 kg
Indicator division:	0.5 kg
Weight:	approx. 3 kg
Dimensions of scales:	approx. 285 x 430 x 83 mm
Distributor: (no service address)	Waagen-Schmitt GmbH Hammer Steindamm 27-29 22089 Hamburg, Germany

Our products are constantly being developed and improved.  
For this reason, design and technical modifications are possible at any time.

## Guarantee

---

The company Waagen-Schmitt GmbH guarantees for 2 years from the purchase date the cost-free rectification of defects on the basis of material or fabrication errors through repair or exchange.

In the event of a guarantee claim, please return the product with the proof of purchase (stipulating the reason for claim) to your dealer.

## Disposal

---

### Disposing of the packaging

Dispose of the packaging with similar materials. Put cardboard with waste paper, films with recyclable materials.

### Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

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